

EasyStand Zing MPS size1



SUPINE TO PRONE
IN ONE MOTION

The **Zing** MPS Advantage

True Multi-Position Stander

- The base allows three positions without changing pads, flipping foot plates or moving the tray
- Provides a smooth transition from a flat-to-load transfer to supine, upright and prone
- Frame will only go supine and upright with Prone Position Lockout engaged

True Hip Pivot Leg Abduction

- Pivots at the hips for proper leg alignment in hips, knees and feet
- Feet maintain angle with the rest of the legs – no wedges needed
- Up to 30 degrees of abduction per leg
- Promotes musculoskeletal development for infants, toddlers and kids

Plenty of Growth

- Size 1 fits children from infancy to 44"
- Can easily grow with the child and adapt to their changing needs
- 70 pound weight capacity



Multi-Use Tray

- Angle adjustable and swings away, the tray opens from either the left or right, or can be removed completely
- Tray can remain attached during transfers
- More activities are possible by swinging tray away during standing therapy



Accommodating Contractures

- No special modifications are needed to accommodate children with contractures
- Using the multi-adjustable foot plates, multi-adjustable knee pads, and the independent calf pads, it is possible to adjust for contractures in one or both sides

Hygienic Covers

- Available for the pelvic, upper body, lateral/hip, head and tray areas
- Easy to wash
- Especially good if stander is used during feeding
- Keeps stander from getting soiled

Unique Lift Mechanisms

- Dual controls are available allowing angle adjustment by foot pedal or hand control
- Gas Spring Lockout prevents an inadvertent change of position, especially important when around other children

Remarkably Modular

- Can be configured as a simple stander or fit a child with the most complex needs
- Designed to grow with the child – options can be added as needed
- Great for a multi-use school or facility



Options



PA5520
Zing MPS Base

Comes with four swivel locking casters, foot operated gas spring, mast receiver with prone position lockout.



PA5608
Directional Locking Caster

Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.



PA5548, PT50080
Foot Straps

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps fit foot height 1.5"-3" (4-8cm)-extra small, 3"-6" (8-15cm)-small.



PA5566, PA5570
Form to Fit Pads

Form to fit pad can be custom fitted to the user. Pad sizes: PA5566-5"H (13cm)-extra small, PA5570-5"H (13cm)-small.



PA3030, PA3032
Hygienic Covers

Washable slip on cover available in six different patterns. For form to fit pads PA5566 and PA5570.



PA5580
Lateral Supports w/Elbow Stop/Arm Rest

Works as a lateral support, elbow rest and armrest. Great for the supine, upright and prone position. For planar only.



PA3039
Hygienic Cover

Washable slip on cover available in six different patterns. For lateral supports PA5576.




PA5596
Black Molded Tray

Tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height. Tray size is 11"Hx19"W (28x48cm).



PA3038
Hygienic Padded Cover

Washable slip on cover is padded and available in six different patterns.

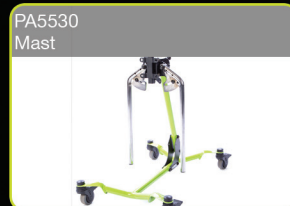

PA5530
Mast

Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head supports.



PA5536
Platform Foot Plate

Adjustable in height. Foot holders and foot straps can be added for more positioning of the feet.



PA5550
Calf Pad with Knee Strap

Gives support from the back of the calf. Adjustable in height and fore/aft.



PA5574
Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm).




PA5562, PA5564, PA5616, PA5618
Planar Pads

Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13-18cm), PA5564-5"Hx9"W (13-23cm), PA5616-7"Hx7"W (18-18cm), PA5618-7"Hx9"W (18-23cm).



PA5630
Elbow Stop with Arm Rest

Great for loading and bringing to upright. Works with both planar and form to fit upper body supports. For planar and form to fit pads.

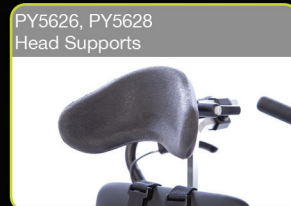


PY5626, PY5628
Head Supports

Extra small head support is 5"Hx8"W (13x20cm) and small head support 6"Hx10"W (15x25cm). Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).



PY5600
Multi-Adjustable Black Molded Tray

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Hx19"W (28x48cm).

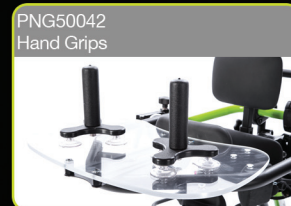

PNG50042
Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups.




PA5532
Mast with Leg Abduction

Mast is the base for positioning options and allows for 0-30° of abduction for each leg.



PA5624, PA5626, PA5628
Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/- .75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot plate sizes: PA5624-6"Lx3"W (15-8cm), PA5626-7.75"Lx3.25"W (20-8cm), PA5628-9.75"Lx4"W (25-10cm).

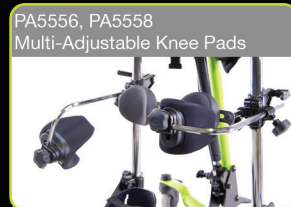

PA5556, PA5558
Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/- .75" (2cm), width 4"-9" (10-23cm), in addition to depth range of 5" (13cm).



PA3011
Positioning Belt

2" (5cm) wide belt has a D-Ring Velcro closure, provides hip stability for user. Fits hip circumference up to 24" (61cm).



PA5568, PA5572, PA5620, PA5622
Form to Fit Pads

Form to fit pad with built in lateral support that can be custom fitted to the user. Pad sizes are: PA5568-5"H (13cm)-extra small, PA5572-5"H (13cm)-small, PA5620-7"H (18cm)-extra small, PA5622-7"H (18cm)-small.



PT50086, PT50088
X-Style Chest Vest

Neoprene® with camlock adjustment provides anterior trunk support for child.



PA5586
Form to Fit Headrest

Head support 6"Hx10"W (15x25cm). Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).



PA5602
Multi-Adjustable Clear Acrylic Tray

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Hx19"W (28x48cm).



PA5594
Push Handle

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child.



PA5612
Dual Control

Allows both foot and hand control of the gas spring lift. Must order push handle PA5594.



PA5540, PA5542, PA5544
Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15-8cm), PA5542-7.75"Lx3.25"W (20-8cm), PA5544-9.75"Lx4"W (25-10cm).



PA5560
Multi-Adjustable Calf Pads

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable up and down with knee pad and fore/aft 4" (10cm) independently.



PA5610
Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits hip circumference up to 24" (61cm).



PA5576
Lateral Supports

Keeps the upper torso in symmetrical alignment. Adjustable from 7-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar only.



PA3028, PA3029, PA3045, PA3046
Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5562, PA5564, PA5616, PA5618.



PA3035, PA3036
Hygienic Covers

Washable slip on cover available in six different patterns. For PY5626 and PY5628.




PA5604
Forearm Wings-Black Molded Tray

Wings give additional area to rest arms when in standing position.



P80809
Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps. Tools not included.

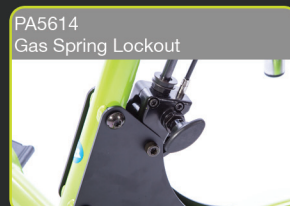

PA5614
Gas Spring Lockout

Locks gas spring in any position from fully supine to prone.



PA5546, PT50082
Secure Ankle/Foot Straps

Velcro® and D-ring adjustment holds feet securely in place at ankle and toes with multiple attachment slots for proper foot positioning. Straps fit foot height 1.5"-3" (4-8cm)-extra small, 3"-6" (8-15cm)-small.



PA5562, PA5564
Planar Pads

Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13x18cm), PA5564-5"Hx9"W (13x23cm).



PA3028, PA3029
Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5562 and PA5564.



PA5578
Lateral Supports with Elbow Stop

Works as a lateral support and elbow rest. Great for the supine position. For planar only.



PA3031, PA3033, PA3047, PA3048
Hygienic Covers

Washable slip on cover available in six different patterns. For form to fit pads PA5568, PA5572, PA5620, PA5622.



PA3037
Hygienic Cover

Washable slip on cover available in six different patterns. For PAC5586.



PA5606
Forearm Wings-Clear Acrylic Tray

Wings give additional area to rest arms when in standing position.



P80246
Angle Locator





Helps therapists monitor standing progress. Attaches to the mast. Can be used on leg abduction.

Health Benefits


- Improve/maintain range of motion
- Decrease joint/muscle contracture
- Management of atrophy in the trunk and leg muscles
- Improve strength to trunk and lower extremities
- Decrease muscle spasms
- Improve/maintain bone integrity/skeletal development
- Lessen/manage the progression of scoliosis
- Manage pressure (ulcers) through changing positions
- Improve bowel function and regularity
- Aid in kidney and bladder functions
- Strengthen cardiovascular system and build endurance
- Improve circulation
- Reduce swelling in lower extremities



Frame Colors

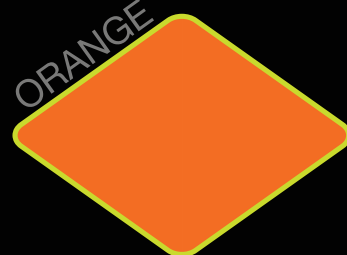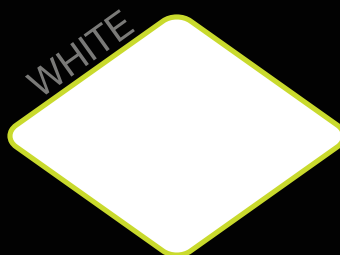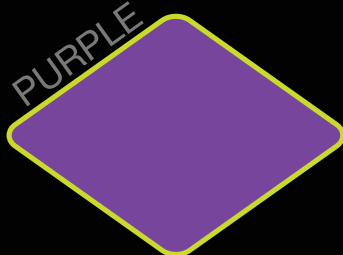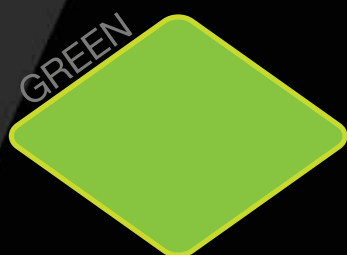

Hygienic Cover Patterns

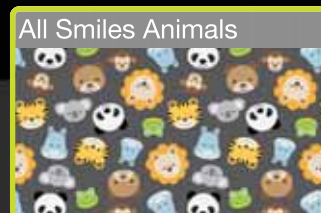

EasyStand

262 West 1st St • Morton, MN 56270, USA • 800.342.8968 • info@easystand.com • easystand.com

Standing technology should only be used under the guidance of a physician with recommendations for standing program protocol and any medical precautions. Standing programs should be monitored by the attending therapist. AMI maintains a policy of continual product improvement and reserves the right to change features, specifications, and prices without prior notification. Check with AMI for latest information. FORM ESZINGBRO-040615 REV A Copyright © 2015 Altimate Medical, Inc. All rights reserved. U.S. & international patents pending. Printed in the U.S.A.